



NEW
PRODUCT

DURABLE, LIGHTWEIGHT
AND PORTABLE



19.8 kg

Total weight without
batteries when folded

Pride[®]
Mobility Products Ltd.
— LIVE YOUR BEST[®] —



Tel: 01869 324600

Email: sales@pride-mobility.co.uk

Web: www.pride-mobility.co.uk

(Version 1: February 2017)


Features and Specification



Introduction

The new i-Go from Pride Mobility Products® Ltd, features advanced folding technology. Enabling it to be quickly transported with ease. The i-Go features front suspension, durable seating system, under seat storage and much more. Designed to fit inside any small space. The i-Go is the perfect choice for the active individual.

Unique Features available with the i-Go

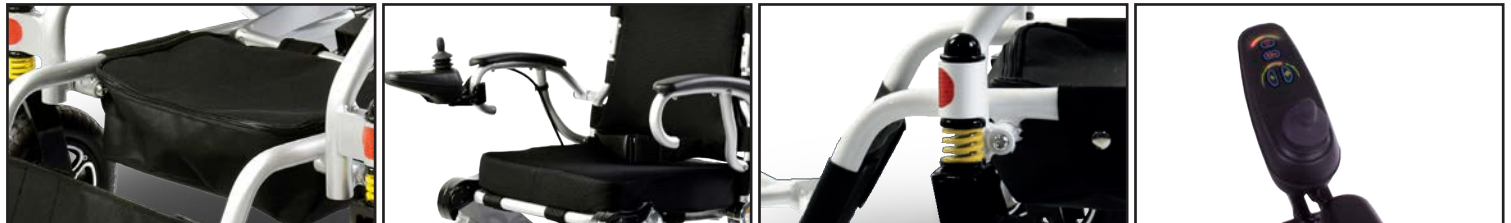

Large secure storage compartment

Unique comfort seating

Hi-torque front suspension

Micon control system



Colours
Only available in Silver

HOW TO FOLD



1 Switch off and remove Joystick



2 Place joystick in backrest pocket and remove cushion.



3 Fold backrest down then pull both arm rests anticlock-wise toward the back.



4 Pull the two black toggles (Located on the back) together.



5 Push down on the seat frame until unit fully flattens.



6 The i-Go is now ready to go...

Technical Specifications

Drive Wheels (cm) 21.06 cm Solid	Caster Wheels (cm) Front 17.08 cm	Max Speed (mph) 3.7 Mph	Ground Clearance (cm) Anti tips 4.4 cm	Turning Radius (cm) 90.00 cm (35.4")	Overall Length (cm) 87.63 cm (34.5")	Overall Width (cm) 47 cm (18.5")	Seating Saddle Back Seat	Batteries x2 6.6 Ah Li-ION battery	Available Electronics 50Ah Micon	Weight Capacity (stone) 18 Stone 8 lbs. (120 k.g)
--	---	---------------------------------------	--	--	--	--	--	--	--	---

Product Weight without Batteries(kg) 19.8 kg	Battery Weight Each (kg) 1.35 kg
--	--

Warranty

- Two-year limited warranty on frame
- One-year limited warranty on Controller
- One-year limited warranty Electronic Components, and Battery Charger.
- 6 month warranty for batteries

Notes

- Speed varies with user weight, terrain type, battery charge, battery condition and tyre pressure.
- Due to manufacturing tolerances and continual product improvement, this specification can be subject to variance of + or - 3%. Dimensions listed are for power base only. Overall measurements will vary based on seating and accessory selections.
- Battery weight may vary based on manufacturer.